




PRODUCT FEATURES:

- Complies with IEC 60336
- Fast and accurate measurements of focal spot size

RaySafe Pro-Slit Camera

The RaySafe Pro-Slit device is a so called slit camera for accurate measurement of the focal spot size according to IEC 60336. Its smart design enables repeatability and accuracy and the possibility to measure the size of any focal spot with one tool.

This test tool can be used with the Pro-Stand (not included) in order to enable an easy and repeatable measurement setup.

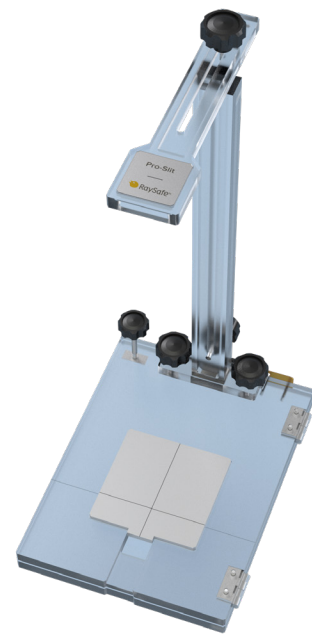

TECHNICAL DATA Product No. 4829300

- Dimensions: 40x40x11.5 mm (1.57x1.57x0.45 in.)
- Material: Tungsten embedded in PMMA
- Slit size: 0.01 mm ($3.9 \cdot 10^{-4}$ in.) wide and 10 mm (0.39 in.) long slit with 8° spread